

CytoGrow™ Recombinant Fetuin A (Bovine)


Catalog Numbers: CG017-10MG, CG017-100MG, CG017-1G

Revision 1.0

Specifications

Synonym	Alpha-2-HS-Glycoprotein; Alpha-2-Z-Globulin; Ba-Alpha-2-Glycoprotein; Fetuin-A; AHSG; FETUA
Species	Bovine
Size	10µg/50µg/1mg
Tag	Tag Free
Purity	≥95%
Endotoxin	≤0.5 EU/mg by the LAL method
Expression System	CHO
Expression Region	I19-I359
pH	7.0-8.0
Appearance	White Powder or clear liquid
Formulation	Lyophilized
Shipping	The product is shipped with blue ice
Reconstitution	It is recommended to redissolve in sterile deionized water.
Storage and Stability	<p>Lyophilized state: 36 months at -20°C to -80°C</p> <p>After reconstitution and under sterile conditions. 6 months at -20°C to -80°C 7-10 days at 2°C to 8°C</p> <p>Store in a manual defrost freezer and avoid repeated freeze-thaw cycles.</p>

Bioactivity

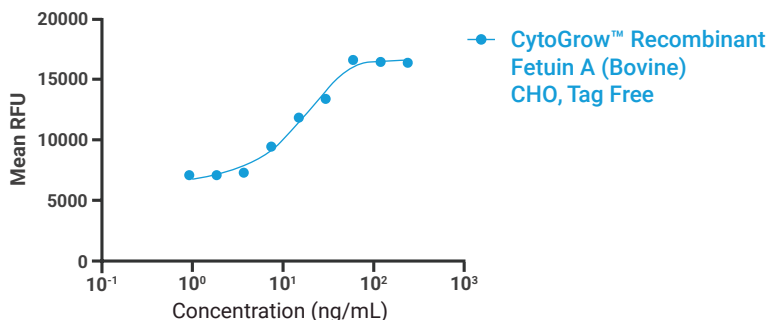

Biological activity is assessed based on the stimulation of B16-F1 cell adhesion. The expected ED₅₀ is less than or equal to 100 µg/mL.

Reconstitution

1. Reconstitute protein to concentrations of 100 µg/mL or higher.
2. Use sterile deionized water to solubilize the lyophilized protein.
3. Gently tap glass vial to help collect the powder to the bottom.
4. Add desired amount of buffer, observe after 2 min. Gentle agitation by hand may be required.
5. Prepare aliquots of solution. Freeze at -20°C to -80°C to avoid repeated freeze thaw cycles.