

CytoGrow™ Recombinant IFN γ (Human)


Catalog Numbers: CG095-A, CG095-B, CG095-C

Revision 1.0

Specifications

Synonym	
Species	Human
Size	10 μ g/vial, 50 μ g/vial, 1 mg/vial
Tag	Tag Free
Purity	$\geq$ 95%, by SDS-PAGE (under reducing (R) & Non-reducing conditions, visualized by Coomassie staining)
Endotoxin	$\leq$ 10 EU/mg by the LAL method
Expression System	E. coli
Expression Region	Gln24-Gln166
pH	7.0-8.0
Appearance	White powder,colorless clear liquid after reconstitution
Formulation	Lyophilized from a 0.22 μ m filtered solution containing 10mM PB,5% Trehalose and 0.01% Tween 80, pH7.4
Shipping	US - Ships overnight
Reconstitution	Recommended to redissolve in sterile deionized water.
Storage and Stability	Lyophilized state: 36 months at -20°C to -80°C After reconstitution and under sterile conditions. 6 months at -20°C to -80°C 7-10 days at 2°C to 8°C Store in a manual defrost freezer and avoid repeated freeze-thaw cycles.

Bioactivity

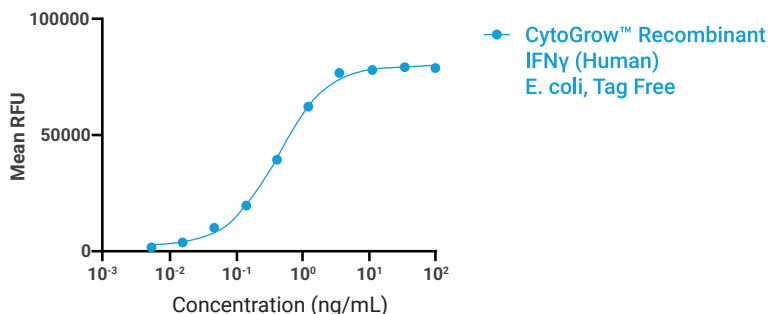

Measured in a cell proliferation assay using Mo7e cells.

Reconstitution

1. Reconstitute protein to concentrations of 100 μ g/mL or higher.
2. Use sterile deionized water to solubilize the lyophilized protein.
3. Gently tap glass vial to help collect the powder to the bottom.
4. Add desired amount of buffer, observe after 2 min. Gentle agitation by hand may be required.
5. Prepare aliquots of solution. Freeze at -20°C to -80°C to avoid repeated freeze thaw cycles.