


CytoGrow™ Recombinant LIF (Bovine)

Catalog Numbers: CG057-A, CG057-B, CG057-C

Revision 1.0

Specifications

Synonym	Leukemia inhibitory factor
Species	Bovine
Size	10µg/50µg/1mg
Tag	Tag Free
Purity	≥95%, by SDS-PAGE (under reducing (R) & Non-reducing conditions, visualized by Coomassie staining)
Endotoxin	≤10 EU/mg
Expression System	E. coli
Expression Region	Ser23-F202
pH	7.0-8.0
Appearance	White Powder or clear liquid
Formulation	Lyophilized
Shipping	US - Ships overnight
Reconstitution	Recommended to redissolve in sterile deionized water.
Storage and Stability	Lyophilized state: 36 months at -20°C to -80°C After reconstitution and under sterile conditions. 6 months at -20°C to -80°C 7-10 days at 2°C to 8°C Store in a manual defrost freezer and avoid repeated freeze-thaw cycles.

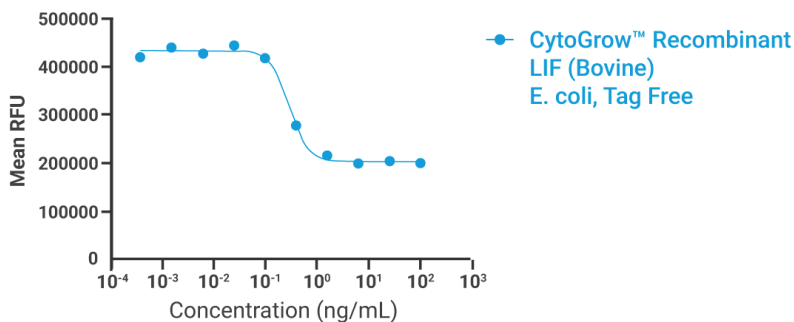

Reconstitution

1. Reconstitute protein to concentrations of 100 µg/mL or higher.
2. Use sterile deionized water to solubilize the lyophilized protein.
3. Gently tap glass vial to help collect the powder to the bottom.
4. Add desired amount of buffer, observe after 2 min. Gentle agitation by hand may be required.
5. Prepare aliquots of solution. Freeze at -20°C to -80°C to avoid repeated freeze thaw cycles.

Biological activity of the CytoGrow™ Recombinant LIF (Bovine), E. coli, Tag Free was measured in the M1 cell differentiation assay.